

DEPARTMENT OF HEALTH CARE ACCESS AND INFORMATION FACILITIES DEVELOPMENT DIVISION

APPLICATION FOR HCAI PREAPPROVAL OF	OFFICE USE ONLY						
MANUFACTURER'S CERTIFICATION (OPM)	APPLICATION #: OPM-0267						
HCAI Preapproval of Manufacturer's Certification (OPM)							
Type: New X Renewal/Update							

Manufacturer Information					
Manufacturer: Getinge IC Production Poland S	Sp. z o.o.				
Manufacturer's Technical Representative: Maciej Przybylak					
Mailing Address: ul. Szkolna 30 62-064 Plewisk	ska Poland, Chino Hills, CA 91709				
Telephone: (909) 606-7622	Email: maciej.przybylak@getinge.com				
NE					
Product Information					

Product Name: LSS-275 Series Steam Sterilizer Vessel Size: 21.5" X 21.5" x 38"	
Product Type: Other electrical and mechanical components	
Product Model Number: LSS-275 and LSS-275 Bio Decon	
General Description: Steam Sterilizer used to sanitize medical instruments, gowns, etc.	

Applicant Information					
Applicant Company Name: EASE LLC.					
Contact Person: Tiffany Tonn					
Mailing Address: 1515 FAIRVIEW AVE, STE 205, MISSOULA, MT 59801					
Telephone: (406) 541-3273	Email: tiffany@easeco.com				
Title: Office Manager					

"Access to Safe, Quality Healthcare Environments that Meet California's Diverse and Dynamic Needs"

14/1/W



STATE OF CALIFORNIA – HEALTH AND HUMAN SERVICES AGENCY



DEPARTMENT OF HEALTH CARE ACCESS AND INFORMATION FACILITIES DEVELOPMENT DIVISION

Registered Design Professonal Preparing Engineering Recommendations					
Company Name: EASE LLC					
Name: Jonathan Roberson California License Number: S4197					
Mailing Address: 5877 Pine Ave., Suite 210, Chino Hills, CA 91709					
Telephone: (951) 295-1892 Email: jon@EASECo.com					
HCAI Special Seismic Certification Preapproval (OSP)					
Special Seismic Certification is preapproved under OSP OSP Number:					
OR CODE O					
Certification Method					
Testing in accordance with: ICC-ES AC156 FM 1950-16					
Other(s) (Please Specify):					
*Use of criteria other than those adopted by the California Building Standards Code, 2022 (CBSC 2022) for component supports and attachments are not permitted. For distribution system, interior partition wall, and suspended ceiling seismic bracings, test criteria other than those adopted in the CBSC 2022 may be used when approved by HCAI prior to testing.					
X Analysis					
Experience Data DATE: 04/05/2023					
Combination of Testing, Analysis, and/or Experience Data (Please Specify):					
COE.					
HCAI Approval					
Date: 4/5/2023					
Name: William Staehlin Title: Senior Structural Engineer					
Condition of Approval (if applicable):					

"Access to Safe, Quality Healthcare Environments that Meet California's Diverse and Dynamic Needs"







5877 Pine Ave. Ste. 210 Chino Hills, CA. 91709 Phn: (909) 606-7622

The Department of Health Care Access and Information PREAPPROVAL OF MANUFACTURER'S CERTIFICATION **OPM-0267**

THIS PREAPPROVAL CONFORMS TO THE 2022 CALIFORNIA BUILDING CODE

GETINGE USA, INC. MANUFACTURER:

Sheet: 1 of 9

EQUIPMENT NAME:

LSS-275 SERIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38" Date: 3/30/23

GENERAL NOTES

- 1. THIS HCAI PREAPPROVAL OF MANUFACTURER'S CERTIFICATION (OPM) IS BASED ON THE 2022 CBC. THE DEMANDS (DESIGN FORCES) FOR USE WITH THIS OPM SHALL BE BASED ON THE 2022 CBC
- 2. THIS DOCUMENT MAY ONLY BE USED WITH THE EXPRESS WRITTEN CONSENT OF THE MANUFACTURER LISTED ABOVE FOR THE SPECIFIC PROJECT SITE AND INSTALLATION LOCATION. THIS DOCUMENT IS INVALID WITHOUT SUCH CONSENT.
- 3. THIS PREAPPROVAL CONFORMS TO THE 2022 CALIFORNIA BUILDING CODE WHERE SDS IS NOT GREATER THAN 1.20, 1.90 & 2.20. SEE DETAIL FOR APPLICABILITY / >
- 4. FORCES PER ASCE 7-16 SECTION 13.3.1, EQUATIONS 13.3-1, 13.3-2 & 13.3-3,
 - WHERE SDS = 1.20, a_D = 1.0, I_D = 1.5, R_D = 1.5, z/b = 0 AT CONCRETE SLAB. SEE FOLLOWING SHEETS FOR Ω_0
 - WHERE SDS = 2.20, \mathbf{a}_{p} = 1.0, \mathbf{I}_{p} = 1.5, \mathbf{R}_{p} = 1.5, \mathbf{z}/h = 0 AT CONCRETE SLAB. SEE FOLLOWING SHEETS FOR Ω_{a}
 - WHERE SDS = 1.90, a_p = 1.0, I_p = 1.5, R_p = 1.5, & z/h < 1 AT CONCRETE SLAB ON METAL DECK. SEE FOLLOWING SHEETS FOR Ω₀
- THIS PREAPPROVAL COVERS ONLY THE SUPPORTS AND ATTACHMENTS OF THE EQUIPMENT TO THE STRUCTURE.
- ALL DESIGN FORCES SHOWN ON THE DRAWINGS ARE FACTORED LOADS THAT SHALL BE USED FOR STRENGTH DESIGN.
- 7. CONCRETE SLAB ON METAL DECK DETAIL VALID FOR DEMANDS SHOWN AT ANY ELEVATION IN THE BUILDING. (i.e. z/h < 1)
- 8. CONCRETE SLAB DETAIL VALID FOR DEMANDS SHOWN AT OR BELOW GRADE. (i.e. z/h = 0)

9. RESPONSIBILITIES OF THE STRUCTURAL ENGINEER OF RECORD OF THE BUILDING

- A. PROVIDE SUPPORTING STRUCTURE TO SUPPORT WEIGHTS AND FORCES SHOWN IN ADDITION TO ALL OTHER LOADS.
- B. VERIFY THAT THE INSTALLATION IS IN CONFORMANCE WITH THE 2022 CBC AND WITH THE DETAILS. MATERIAL AND GAGE OF THE UNIT WHERE ATTACHMENTS ARE MADE AGREE WITH THE INFORMATION SHOWN ON THE PREAPPROVAL DOCUMENTS.
- C. VERIFY THAT PROJECT SPECIFIC VALUES OF SDS & Z/h RESULT IN SEISMIC FORCES (Eh, Ev) THAT DO NOT EXCEED THE VALUES ON THE DETAILS.
- D. VERIFY THAT THE CONCRETE SLAB TO WHICH THE EQUIPMENT IS ANCHORED MEETS THI REQUIREMENTS OF THE APPLICABLE ICC ESR REPORT, AND THIS OPM.
- E. VERIFY THAT THE ANCHORS ARE AN ADEQUATE DISTANCE FROM ANY SLAB EDGES OR OPENINGS (SEE TYPICAL DETAIL ON SHEET 2).
- F. VERIFY THAT ALL NEW OR EXISTING ANCHORS ARE AN ADEQUATE DISTANCE FROM THE UNIT ATTACHMENTS AND CHECK FOR INTERACTION WHERE OTHER ANCHORS ARE WITHIN 18" OR 6hef FROM THIS UNIT'S ANCHORS.



www.EquipmentAnchorage.com

OF

GETINGE USA, INC.

DES. J. ROBERSON

JOB NO. 36-2301

DATE

3/30/23

SHEET 2

9 SHEETS

LSS-275 SÉRIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38"


10. EXPANSION ANCHORS:

A. ATTACHMENT IS TO BE MADE WITH THE ANCHORS LISTED BELOW AND INSTALLED AS DESCRIBED IN THE CORRESPONDING ICC REPORT.


Anchor Diameter	Concrete Type	Min. f'c (psi)	Anchor Type	ICC Report No.	Min. Embed.	Min. Spacing	Min. Edge Dist.	Min. Conc. Thickness	Torque Test	Direct Tension Test
3/8"	Sand Light Weight	3000	Hilti Kwik Bolt TZ2	ESR-4266	2"	6.75"	12"	See Detail "A"	30 FT-LB	N/A
1/2"	Normal Weight	3000	Hilti Kwik Bolt TZ2	ESR-4266	2"	6"	14"	4"	40 FT-LB	1605 lb
1/2"	Normal Weight	3000	Hilti Kwik Bolt TZ2	ESR-4266	3.25"	6"	24"	6"	40 FT-LB	2685 lb

- B. THIS PREAPPROVAL ALLOWS FOR UP TO A MAXIMUM OF 2 ADJACENT CONCRETE SLAB EDGES, SEE SCHEDULE AWAY MINIMUM (i.e. CORNER).

 SEE ADJACENT DETAIL FOR ADDITIONAL MINIMUM ALLOWABLE CONCRETE EDGE DISTANCES.
- C. TESTING AND SPECIAL INSPECTION OF EXPANSION ANCHORS SHALL
 BE PERFORMED BY AN APPROVED INDEPENDENT AGENCY
 EMPLOYED BY THE FACILITY OWNER PER CBC 1704A & 1910A.5
 AND CAC 7-149. ALL REPORTS SHALL BE SENT TO THE INSPECTOR
 OF RECORD, OWNER AND THE ARCHITECT OR ENGINEER IN STARBLE RESPONSIBLE CHARGE.
 - (i) AFTER AT LEAST 24 HOURS HAVE ELAPSED SINCE INSTALLATION,
 DIRECT PULL TENSION TEST OR TORQUE TEST AT LEAST 50% OF
 THE ANCHORS.
 - (ii) ACCEPTANCE CRITERIA:
 - DIRECT TENSION TEST: THE ANCHOR SHOULD HAVE NO
 OBSERVABLE MOVEMENT AT THE TEST LOAD. A PRACTICAL WAY
 TO DETERMINE OBSERVABLE MOVEMENT IS THAT THE WASHER
 BECOMES LOOSE.
 - TORQUE TEST: THE APPLICABLE TORQUE MUST BE ACHIEVED WITHIN THE FOLLOWING LIMITS: WEDGE TYPE: 1/2 TURN OF THE NUT
 - (iii) IF ANY ANCHOR FAILS, TEST ALL ANCHORS.
- D. AVOID DAMAGING EXISTING STEEL REINFORCING IN CONCRETE SLAB WHEN INSTALLING CONCRETE EXPANSION ANCHORS.
- E. PROVIDE FOR FULL THREAD ENGAGEMENT OF NUT & WASHER.
- 11. BOLTS THROUGH CONCRETE ON METAL DECK
 - A. BOLTS SHALL BE TORQUED BY 3/4 TURN OF THE NUTS AFTER THE SNUG TIGHT (THE SNUG-TIGHT CONDITION IS DEFINED AS THE TIGHTNESS REQUIRED TO BRING THE CONNECTED PLIES INTO FIRM CONTACT) CONDITION IS ACHIEVED, UNLESS OTHERWISE NOTED.
 - B. THROUGH BOLT HOLES SHALL BE 1/16" LARGER THAN BOLT SIZE (HOLE SIZE = BOLT SIZE + 1/16) FOR CONCRETE.
 - C. THROUGH-BOLTS IN CONCRETE SHALL RECEIVE SPECIAL INSPECTION AND TESTING (THROUGH BOLTS WITH STEEL TO STEEL CONNECTION IN TENSION DO NOT REQUIRE TENSION TESTING) IN ACCORDANCE WITH REQUIREMENTS FOR POST-INSTALLED ANCHORS.

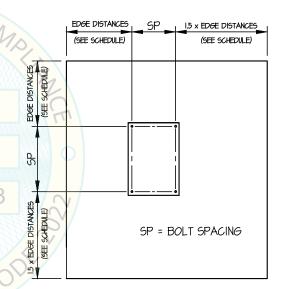

TYPICAL CONCRETE EDGE DETAIL (SLAB ON GRADE ONLY)

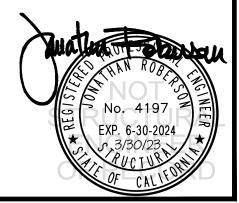

DES. J. ROBERSON

www.EquipmentAnchorage.com

GETINGE USA, INC.


36-2301 JOB NO.

SHEET

3/30/23 DATE

SHEETS

SS-275 SÉRIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38"

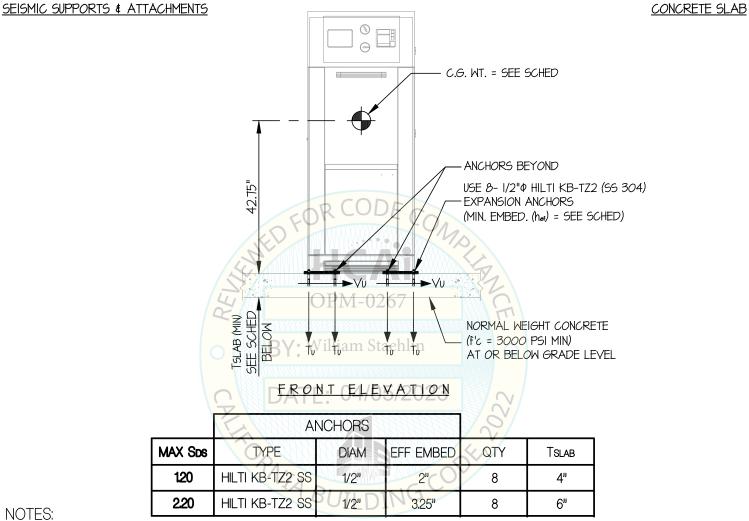

- 1. FORCES ARE DETERMINED PER 2022 CALIFORNIA BUILDING CODE AND ASCE 7-16. STRENGTH DESIGN IS USED. (EXAMPLE: $\Delta p = 1.0$, Ip = 1.5, Rp = 1.5, $\Omega_0 = 2.0$, z/h = 0)
- 2. THIS CALCULATION ENCOMPASSES WEIGHTS AND VERTICAL C.G. POSITIONS NOT EXCEEDING VALUES SHOWN.
- 3. THIS CALCULATION WAS PREPARED WITHOUT KNOWLEDGE OF ANY SITE CONDITION, COMPATIBILITY FOR USE WITH A SITE SHALL BE EVALUATED BY THE STRUCTURAL ENGINEER OF RECORD OF THE INSTALLATION (SEOR). USE REQUIRES APPROVAL BY THE SEOR.
- 4. STRUCTURAL ENGINEER OF RECORD FOR THE INSTALLATION SHALL VERIFY ALL CONDITIONS, EVALUATE INTERACTION WITH ADJACENT EQUIPMENT AND ANCHORS. AND PROVIDE SUPPORT STRUCTURE DESIGNED TO SUPPORT WEIGHTS AND FORCES SHOWN IN COMBINATION WITH ALL OTHER LOADS THAT MAY BE PRESENT.
- 5. SEE GENERAL NOTES: SHEETS 1 AND 2.



www.EquipmentAnchorage.com

GETINGE USA, INC.

SS-275 SÉRIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38"


DES. J. ROBERSON


36-2301 JOB NO.

3/30/23 DATE

SHEET

OF SHEETS

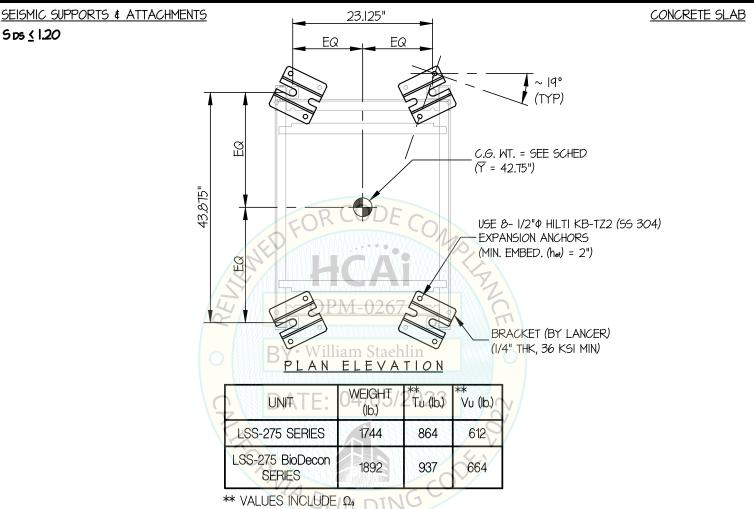

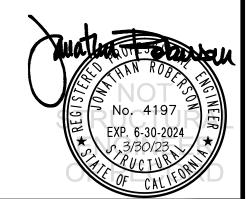

S DS ≤ 1.20

www.EquipmentAnchorage.com

GETINGE USA, INC.

SS-275 SÉRIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38"


DES. J. ROBERSON

36-2301 JOB NO.

3/30/23 DATE

SHEET

SHEETS

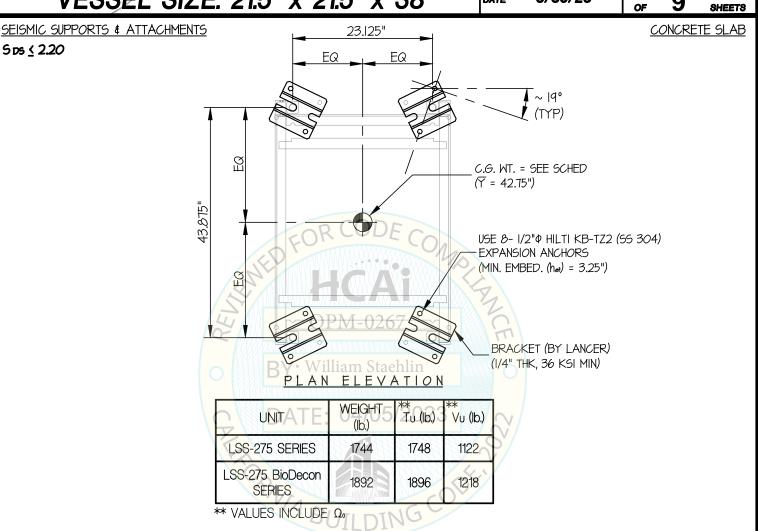



S DS ≤ 2.20

www.EquipmentAnchorage.com

GETINGE USA, INC.

SS-275 SÉRIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38"

DES. J. ROBERSON

JOB NO. 36-2301


DATE 3/30/23

SHEET 6

9 SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

CONCRETE SLAB ON METAL DECK

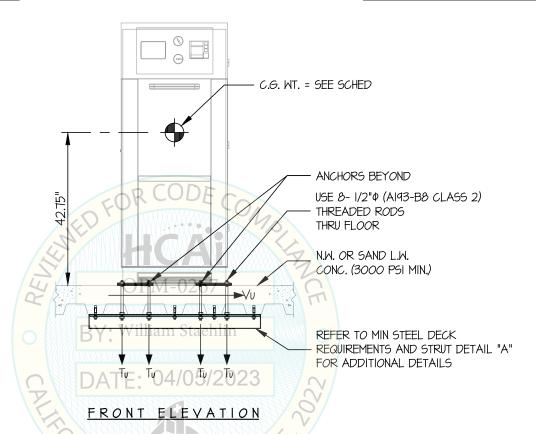

NOTES:

1. FORCES ARE DETERMINED PER 2022 CALIFORNIA BUILDING CODE AND ASCE 7-16. STRENGTH DESIGN IS USED. (EXAMPLE: SDS = 1.90, ap = 1.0, lp = 1.5, Rp = 1.5, Ω_0 = 2.0, z/h \leq 1)

HORIZONTAL FORCE (En) = 2.28 Wp

HORIZONTAL FORCE (Emh) = 4.56 Wp (FOR CONCRETE ANCHORAGE)

VERTICAL FORCE (Ev) = 0.38 Wp

- 2. THIS CALCULATION ENCOMPASSES WEIGHTS AND VERTICAL C.G. POSITIONS NOT EXCEEDING VALUES SHOWN.
- 3. THIS CALCULATION WAS PREPARED WITHOUT KNOWLEDGE OF ANY SITE CONDITION. COMPATIBILITY FOR USE WITH A SITE SHALL BE EVALUATED BY THE STRUCTURAL ENGINEER OF RECORD OF THE INSTALLATION (SEOR). USE REQUIRES APPROVAL BY THE SEOR.
- 4. STRUCTURAL ENGINEER OF RECORD FOR THE INSTALLATION SHALL VERIFY ALL CONDITIONS, EVALUATE INTERACTION WITH ADJACENT EQUIPMENT AND ANCHORS, AND PROVIDE SUPPORT STRUCTURE DESIGNED TO SUPPORT WEIGHTS AND FORCES SHOWN IN COMBINATION WITH ALL OTHER LOADS THAT MAY BE PRESENT.
- 5. SEE GENERAL NOTES: SHEETS 1 AND 2.



www.EquipmentAnchorage.com

OF

GETINGE USA, INC.


SS-275 SERIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38" DES. J. ROBERSON


JOB NO. 36-2301

DATE 3/30/23

SHEET 7

9 SHEETS

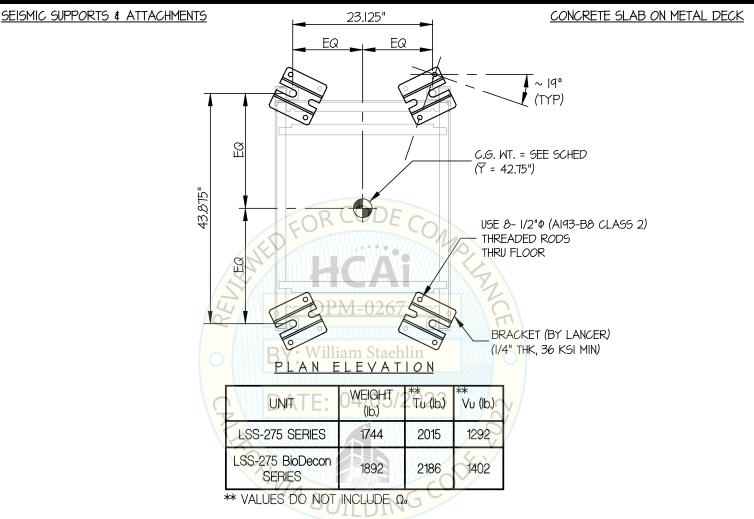

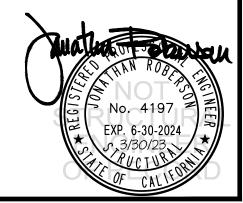

www.EquipmentAnchorage.com

OF

GETINGE USA, INC.

LSS-275 SÉRIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38"


DES. J. ROBERSON

JOB NO. 36-2301

DATE 3/30/23

8HEET

9 SHEETS

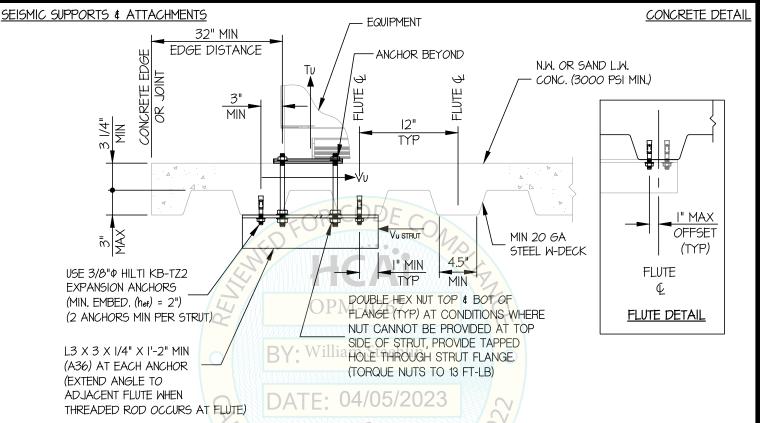

DECK REQUIREMENTS AND STRUT DETAIL

PNIA BUILDING COD



www.EquipmentAnchorage.com

OF

GETINGE USA, INC.

SS-275 SERIES STEAM STERILIZER VESSEL SIZE: 21.5" x 21.5" x 38" DE8. J. ROBERSON

JOB NO. 36-2301

DATE 3/30/23

SHEET

9 знеетз

